

# Sensitive adductor tendon

Example 1



Example 2



Example 3



# With Your® Liner

Limb



3D model



Your® Liner




## Indication

Patient complains about the compression applied by the liner on his medial adductor tendon



## Typical issues

**Off-the-shelf liner applies excessive pressure on adductor tendon**  
Off-the-shelf femoral liners are by essence cylindrical or conical. As off-the-shelf liners are designed to fit residual limbs of various lengths, they do not include a radius at the level of the adductor tendon, and thus apply extra pressure in this region, which can lead to discomfort, or skin irritations.

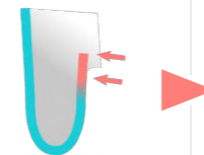

## Key features



**Anatomically tailored**



**Adapted compressions**



**Your® Liner includes radius at adductor tendon to relieve pressure**

Your® Liner being designed for individual patients, a radius can be included in the 3D model to ensure the pressure applied on the adductor tendon is relieved, thus reducing the risk of the liner generating discomfort or skin irritations in this area.

**How can Your® Liner help ?**

